

Q-Switched Nd:YAG

Enables Advanced Combination Treatments



7 spot sizes for customized treatments, 4 of which are available also as disposable tips



2mm




2.5mm



3.5mm




5mm

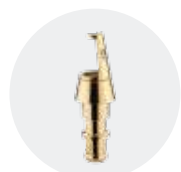


2mm




2.5mm

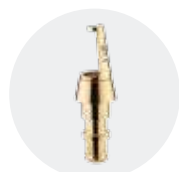

3.5mm

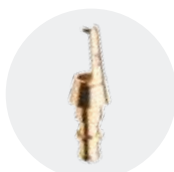


4mm

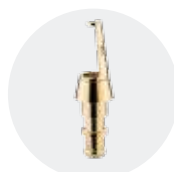

5mm




6mm

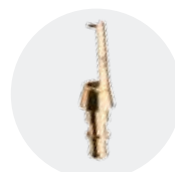


8mm

“Top-Hat” beam profile advantages

- Enhanced treatment safety due to avoidance of bulk heating using the homogenizer
- Top-hat beam profile minimizes epidermal damage, tissue textural changes and the occurrence of scarring
- A lack of an homogenizer might result in hot spots within the beam and increased risk of scarring

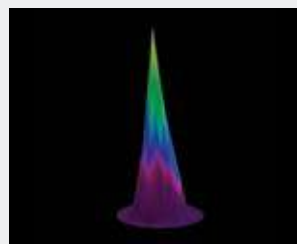

Beam profile side view without an homogenizer

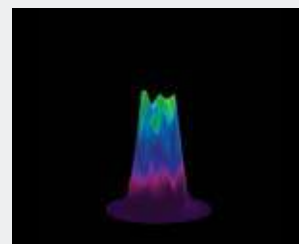

Beam profile side view with an homogenizer for a top-hat beam

For skin toning by treatment of pigmented lesions, dark tattoo removal and advanced combination treatments with other M22 modules

Typical Clinical Results



Courtesy of Hiroki Kanda, MD.



Courtesy of Yuanhong Li, MD.



Courtesy of Michael Oeser, MD.



Courtesy of Moshe Lapidoth, MD.


Specifications

Q-Switched Nd:YAG	
Wavelength	1064 nm
Fluence	0.9-14 J/cm ²
Pulse Duration	6-8 ns

Repetition Rate	0.5-5.0 Hz
Spot Size	2, 2.5, 3.5, 4, 5, 6 and 8 mm

Warnings and risks

The M22™ Q-Switched Nd:YAG is contraindicated for patients with skin disorders or inflammatory skin conditions. The use of the Q-switch could cause redness, swelling, pigmentation changes. See the system user manual for a complete list of contraindications and risks.

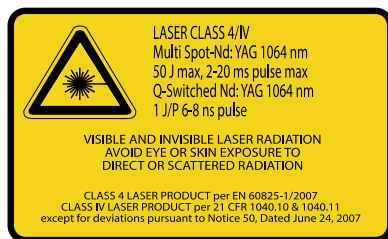

Lumenis Ltd.
Yokneam Industrial Park
6 Hakidma Street P.O.B. 240
Yokneam 2069204, Israel
Tel: +972-4-959-9000

Lumenis (Germany) GmbH
Heinrich-Hertz-Str 3 D-63303
Dreieich-Dreieichenhain
GERMANY
Tel: +49 (0) 6103 8335 0

EMEA
Dreieich Dreieichenhain,
Germany
T +49 6103 8335 0
F +49 6103 8335 300
Roma (RM), Italy
T +39 06 90 75 230
F +39 06 90 75 269
Hertfordshire, UK
T +44 20 8736 4110
F +44 20 8736 4119

AMERICAS
San Jose, CA, USA
T +1 408 764 3000
+1 877 586 3647
F +1 408 764 3999

JAPAN
Tokyo, Japan
T +81 3 6743 8300
F +81 3 6743 8301

ASIA / PACIFIC
Beijing, China
T + 86 10 5737 6677
Gurgaon, India
T + 91 124 422 07 95
Kowloon, Hong Kong
T + 852 217 428 00
F + 852 272 251 51



Rx ONLY

© 2017 All Rights Reserved. The Lumenis Group of Companies. PB-1003630 Rev D. Lumenis, its logo, M22, are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice Lumenis is ISO 13485:2012 certified. | M22 is CE approved. | M22 is FDA cleared.

For more information: www.lumenis.com



LUMENIS.COM/AESTHETIC